

ÖSSUR FORMFIT® PRO ANKLE

Össur Formfit® Pro Ankle is a compressive ankle sleeve with 2 integrated malleoli pads, designed to be worn after mild ankle injuries. The dynamic 3-D MotionTech™ knit is supportive and meets certified medical grade compression requirements. Proven superior to leading brands for moisture wicking and breathability, your ankle stays cool during wear. CoolVent™ technology in the area in front of the ankle further enhances breathability and comfort.



MotionTech™



CoolVent™



SmartPad™



FoldLock™



INDICATIONS FOR USE

- For ankle conditions that may benefit from compression, such as:
- Ankle injury rehabilitation
- Mild sprains and strains
- Mild instability


CONTRAINDICATIONS

- Must not be used by individuals for whom compression is contraindicated.

FEATURES AND BENEFITS

- MotionTech™ for compression and dynamic fit with superior moisture wicking and breathability, keeping your ankle cool during wear
- Certified medical grade compression (20-30 mmHg)
- CoolVent™ technology in the area in front of the ankle for added breathability and comfort
- 2 integrated malleoli SmartPads™ for targeted support and cushioning
- Patented FoldLock™ system for easier and more secure application

REQUIRED MEASUREMENTS

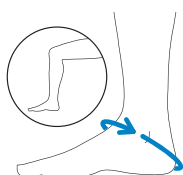

Measure ankle circumference diagonally over the base of the heel and across the ankle bone. Ankle in 90 degrees.

Size	Ankle circumference cm
XS	18 - 20
S	20 - 23
M	23 - 26
L	26 - 29
XL	29 - 32

ÖSSUR FORMFIT® PRO ANKLE

ORDERING INFORMATION

BLUE

Size	Left	Right	Size	Left	Right
XS	FP201201	FP201101	XS	FP200201	FP200101
S	FP201202	FP201102	S	FP200202	FP200102
M	FP201203	FP201103	M	FP200203	FP200103
L	FP201204	FP201104	L	FP200204	FP200104
XL	FP201205	FP201105	XL	FP200205	FP200105

BLACK

